


VISU Pro Standard Model 3

Worktop ceramic, optional compressed air tap, polypropylene interior, vertical sash

VISU Pro Fume Cupboard enables the safe processing of harmful materials and gases.

- **High quality** - Materials for professional use and finished design, almost maintenance-free.
- **Ergonomic** - Clear and optimal design of details.
- **Space saving** - The largest interior volume of the collection in relation to the external dimensions.
- **Safe** - The air suction is on the full width of the fume cupboard on upper and lower edge of the inner back wall.
- **Ecological** - Designed and manufactured responsibly in Finland, made to last for decades.

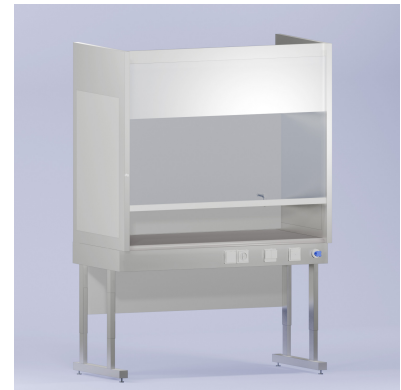

The fume cupboard can be equipped with a compressed air tap. Details are provided below.

Product information

Standard widths	1200, 1500 and 1800 mm
Worktop depth	600 and 800 mm
Total depth	730 and 930 mm
Total height	2680 mm when 1-part sash totally opened
Worktop height adjustment	Manual 750-900 mm, set to 900 mm
Worktop material	Ceramic
Side walls and back panel material	Polypropylene
Front sash	Vertical laminated glass
Frame	Powder coated steel
Ductwork	High-pressure laminate + back panel material
Standards	According to SFS-EN 14175-2 : 2003 safety features: lift limiter in the front sash, set to 500 mm. Sash emergency brake.

VISU KALUSTE OY

📍 Eeliksentie 6, 34870 Visuvesi, Finland

✉️ visu@visu.fi

☎️ +358 44 749 4940

🌐 visu.fi



INNOSTAVIA YMPÄRISTÖJÄ
TIETEELLE JA TUTKIMUKSELLE

Sockets	IP44, 2-socket outlet, 1 pcs
Air flow control	0-1 control switch for ventilation (building services)
Tap (optional)	Compressed air (CA) tap, control valve in front panel.
Residual current device	Included

Additional technical information

The technical properties depend on the fume cupboard width.

Widths, [mm]	1200	1500	1800	
Internal work area widths, [mm]	1130	1430	1730	
Air volume, [l/s]	170	215	260	Sash opened 300 mm, 0.5 m/s inflow, pressure drop 100 Pa.
Standard lights, led W	18	24	30	LUX: 1022, CRI: 80, frequency: 3.3 kHz, color temperature: 4000K
Ventilation connection, [mm]	200	250	2×200 (cc 900 mm)	

Brochures and files

Selection of a suitable fume cupboard

See the table below to select different properties. We always equip the fume cupboards according to customers' needs. Feel free to contact us, and we will design suitable solutions together!

Delivery and installation

Fume cupboards are delivered, if so agreed, installed in place, ready to be connected to building services systems. A fume cupboard must be connected to a separate exhaust air duct. The water supply, drainage, gas and electricity work must be performed by professionals.

VISU KALUSTE OY

Eeliksentie 6, 34870 Visuvesi, Finland

visu@visu.fi

+358 44 749 4940

visu.fi



INNOSTAVIA YMPÄRISTÖJÄ
TIETEELLE JA TUTKIMUKSELLE



VISU KALUSTE OY

📍 Eeliksentie 6, 34870 Visuvesi, Finland

✉️ visu@visu.fi

☎️ +358 44 749 4940

🌐 visu.fi